


Sanipull Assembly Instructions

Sanipull Quick Installation Guide for Bedside Nurse Call Pulls

Scissors are needed to cut the Sanipull ribbon to optimal length.


Attaching the Sanipull ribbon to the adapter

Withdraw from the bulk spool or coil sufficient Sanipull ribbon and cut to the optimum length. From the direction of the middle of the adapter insert the ribbon end under the plastic bridge found in the lower clamshell (look for guide arrows). Fold the ribbon back toward the adapter's bottom end. Close the lower clamshell door.

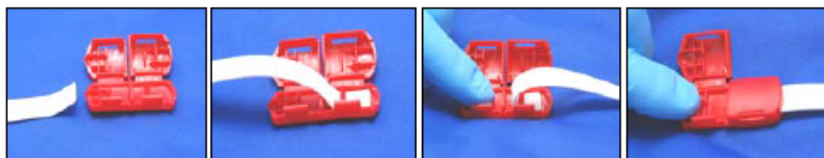

Installation of Sanipull onto nurse call fixture/light pull


1. Cut the nurse call fixture's original pull string/cord about 4" below the wall fixture. Tie a Figure 8 knot at the end of it. Cut off excess. To resist pull through, a #6 stainless steel flat washer can be threaded onto the pull string before making the knot. Special hardware available at no cost.
2. Position the pull cord knot into the top clamshell, ensuring that it sits below the V-shaped notch. Close the upper clamshell door to secure.

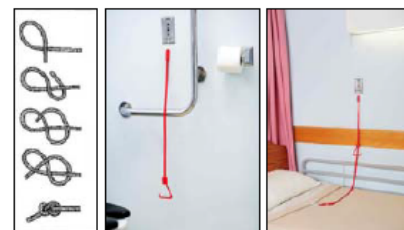

Note: for Austco™ pull cords the fixture's plastic snap connector must be cut off the end of the clear plastic lead. Melt the end of the lead to make a ball end. Insert the ball end into the top end of the clamshell adaptor, taking care that it sits below the V-shaped notch. Close the upper clasp door.

Attaching a Sanipull Bangle onto the Sanipull ribbon


1. Insert the end of the ribbon into the slot at the top of the handle cover (install 2 covers, if 2 handles are used).
2. With the end of the ribbon positioned flush into the ribbon-profile recess of the bottom prong, pass ribbon over the next prong and weave under and over the next two prongs. (Ridges on one side of the prongs allow only one weave direction). A second bangle can be installed mid-ribbon to meet current legislative guidelines.
3. Slide the cover down over the ribbon-threaded posts to secure the grab handle.

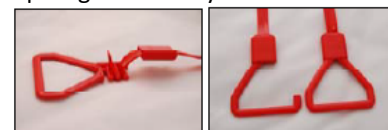

There are two styles of grab handle available for optimum accessibility: open (Bangle- 0,R,W) and closed (BangleC- 0,R,W) Handles with covers are packaged in bags of 10 sets. (0,R and W indicate colour Orange, Red or White)

Attaching a Super Sanipull spring clip to user-end of a Sanipull ribbon

1. Hold the clip upside down. Position the end of the ribbon over the centre of the clip pointing toward the clasp. Following the guide arrow, insert the ribbon under the bridge until the end butts against the end of the clasp.
2. Fold the ribbon back toward the clasp end and close the clamshell door.
3. Attach spring clip onto linen, clothing, around bed rail or wheelchair arm rest.

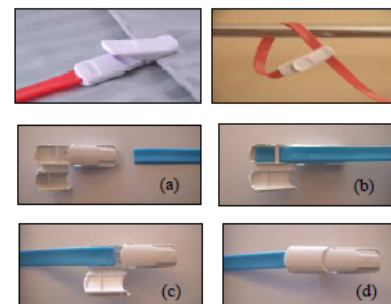

Sanipull is wipeable, to clean use non-caustic alcohol-based sanitiser/gel. Sanitize Sanipulls after installation and regularly according to facility housekeeping policy.